

amalgatome

MD



# Skin Recovery Device Training

*FOR CADAVERIC USE ONLY*

amalgatome | MD

# Device Components



# Transportation Case

Keep this gray box, it will be used to return the unit after trial or for service.



# Complete Amalgatome MD System




# Power Supply & Cord


Does not  
go into  
autoclave.



# Sterilization Tray

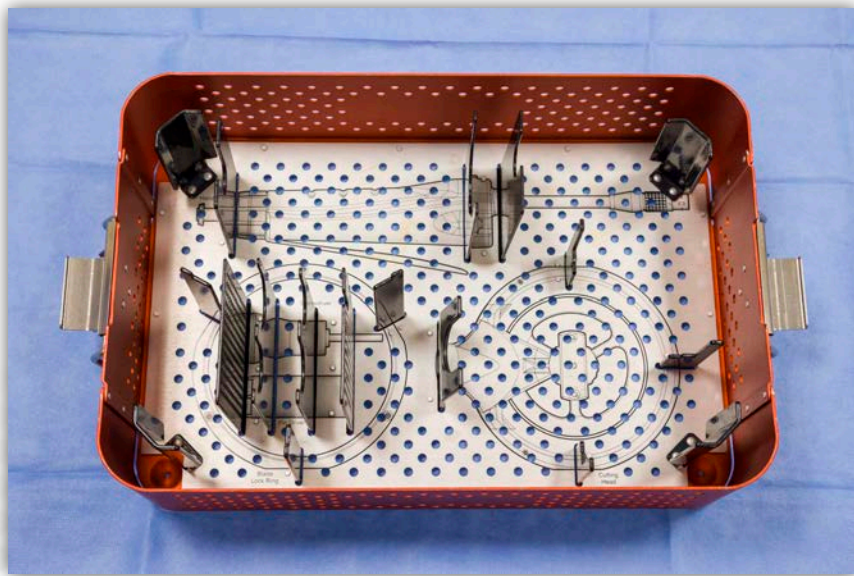


# Sterilization Tray

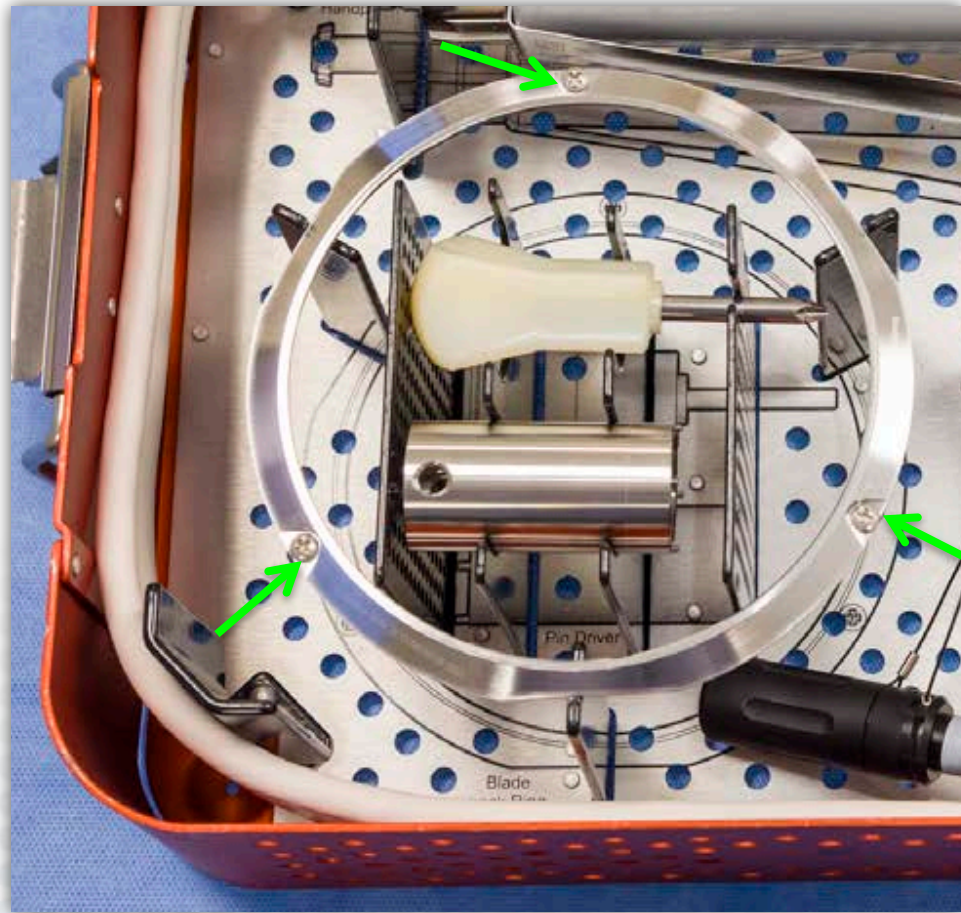

Ensure the blade locking ring is placed inside the tray with the screw heads facing upward.




# Sterilization Tray

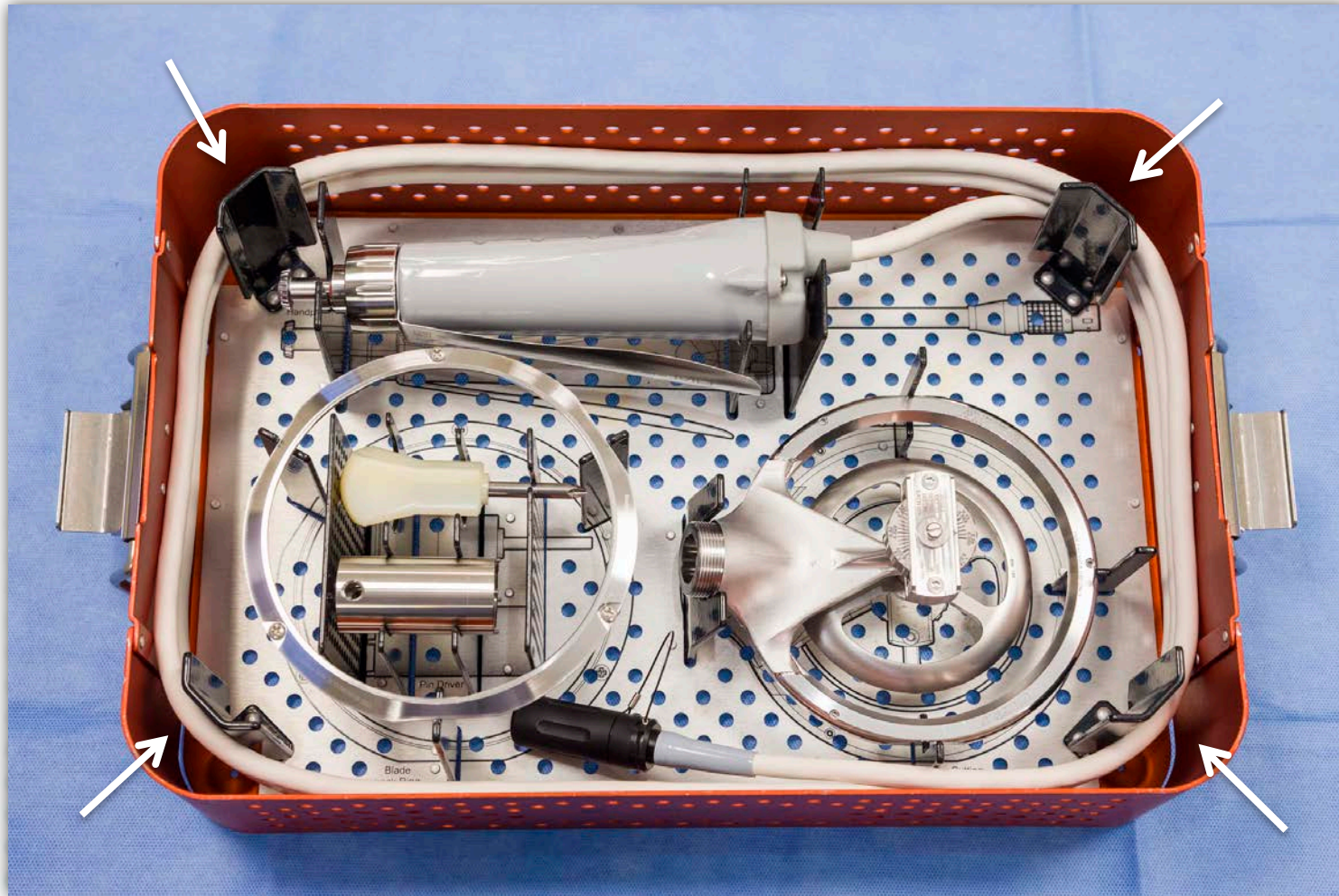

Note: as the power cord is wrapped around the outside of the tray, make sure not to apply too much tension at the corners of the tray.



# Head Assembly

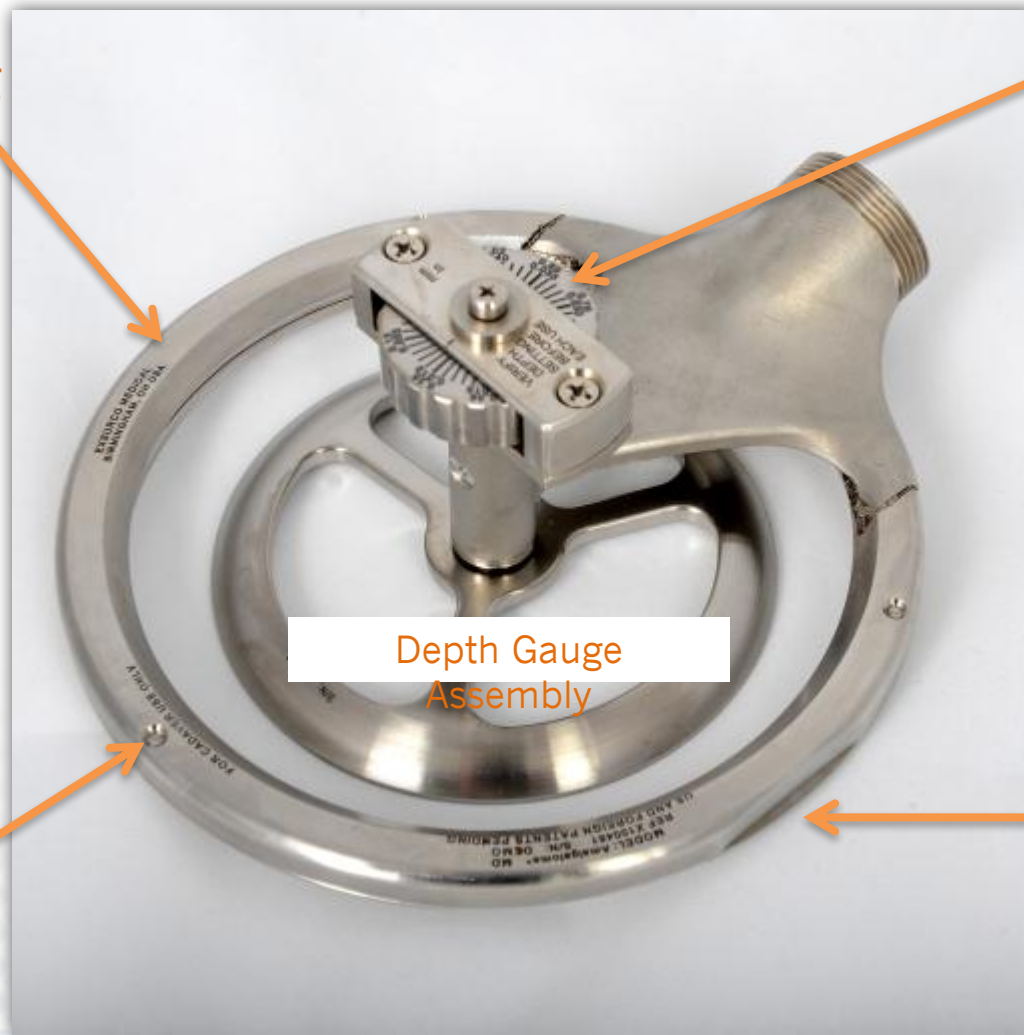

Housing  
Frame

Depth  
Adjustment  
Knob

The head  
assembly is  
not to be  
disassembled.

Depth Gauge  
Assembly

Screw  
Hole

Blade Locking  
Ring w/ 3  
Captive Screws



# Hand Piece

Power Cord  
with  
Protective  
Cap

Activation  
Button


Activation  
Lever

Locking  
Knob

Pinion

Locking Nut

Handle

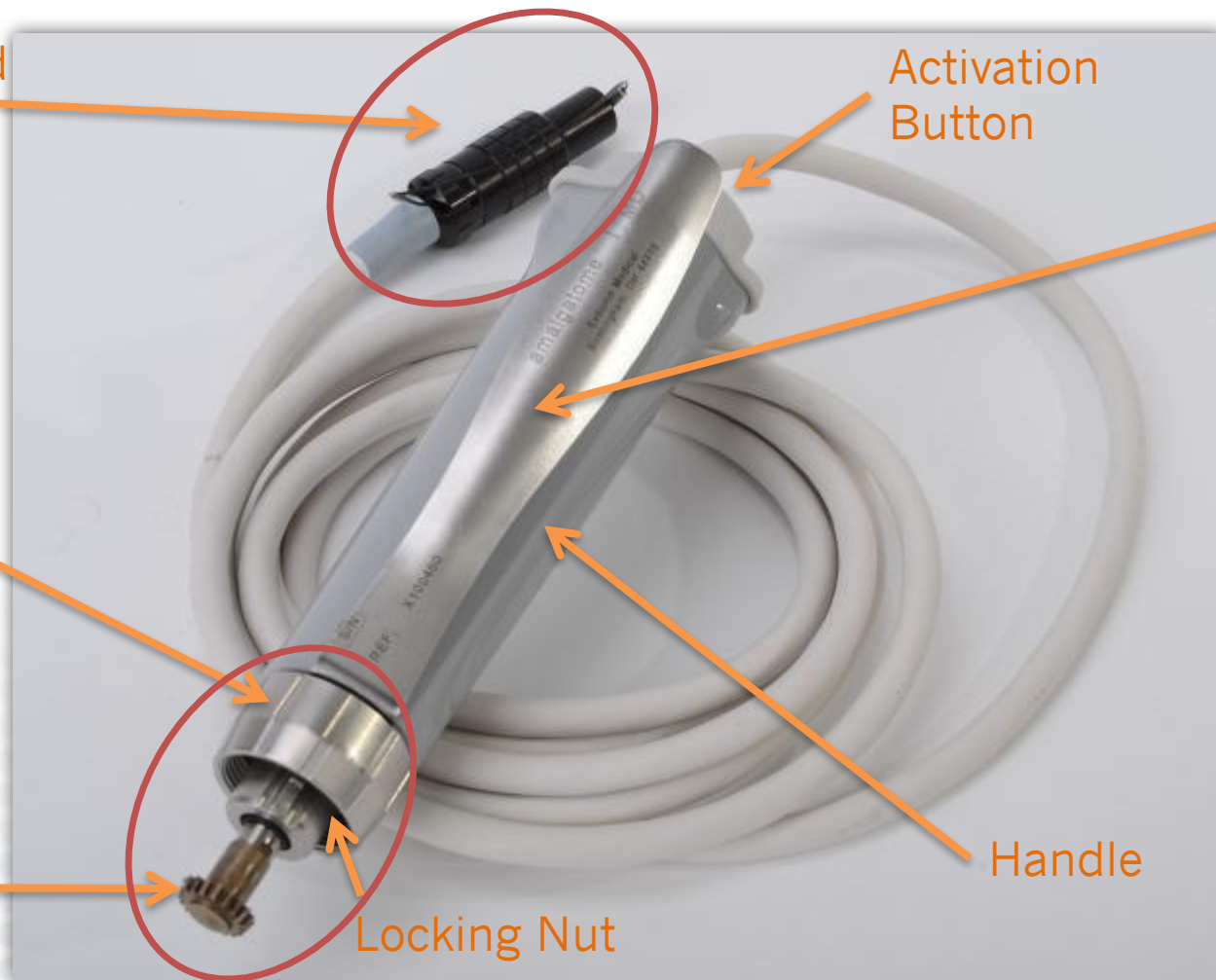

# Disposable Blades

Note:  
Always  
check the  
sterile blade  
package  
integrity  
and  
expiration  
dates  
before use.



Sterile,  
Single  
Use

# Disposable Blades

Please note: this is best practice for handling the blade.



# Disposable Blades

Please note: it is not safe to handle the blade anywhere on the inside of the blade as shown here.



# Operator & Circulator

## Operator



When the term "sterile operator" is used this refers to the technician who has followed protocols for scrubbing, gowning and gloving and is practicing sterile technique.

## Circulator



When the term "circulator" is used this refers to the technician who is not scrubbed in and is providing support within the recovery room.

# Need to get a hold of us?

**Phone: 1-800-243-6049**

Press the number below for the following departments:

- 1 – Sales & Clinical Support
- 2 – Marketing & Product Inquiries
- 4 – Customer Service or to Report a Quality Concern
- 5 – Technical Services
- 6 – Quality Department
- 7 – Accounting Department



**Website:** [www.exsurco.com](http://www.exsurco.com)

**Emails:**

- Sales: [sales@exsurco.com](mailto:sales@exsurco.com)
- Information:  
[info@exsurco.com](mailto:info@exsurco.com)
- Support:  
[support@exsurco.com](mailto:support@exsurco.com)

